



بِسْمِ اللَّهِ الرَّحْمَنِ الرَّحِيمِ

# Biotechnology

## Section 3

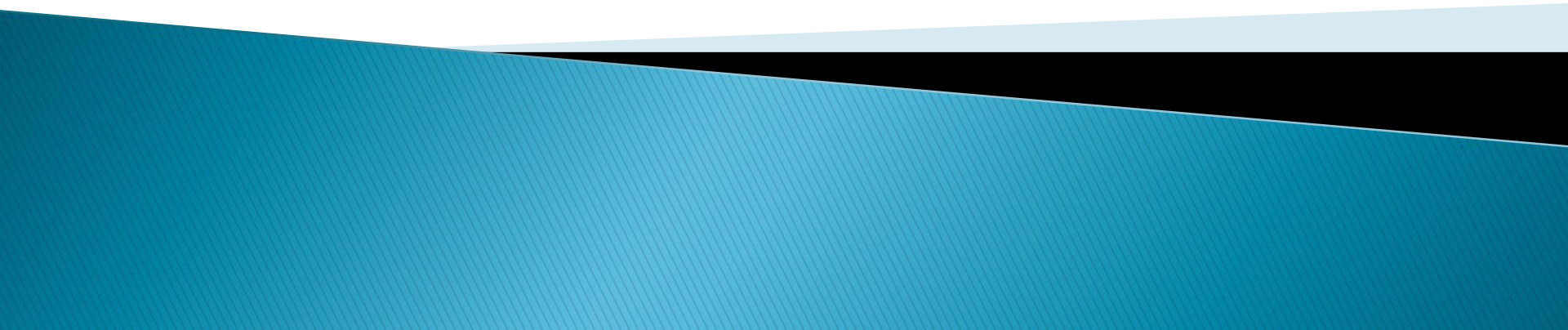

# Objectives

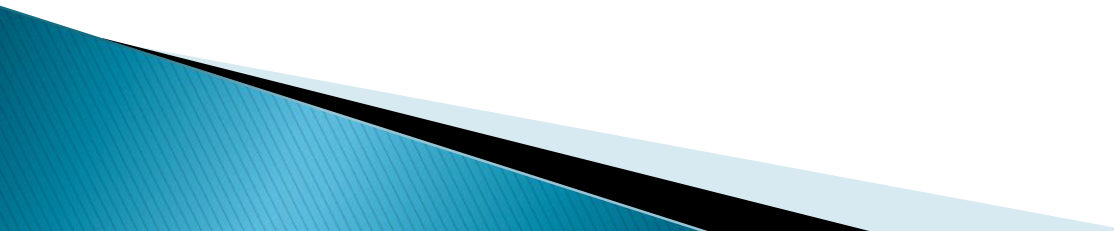- ❑ Macroscopic examination of *Streptomyces*.
  - ❑ Microscopic examination.
  - ❑ Streaking for isolation.
- 

# Macroscopic and Microscopic characters of Streptomyces

- ▶ *Streptomyces* are Gram positive, spore-forming bacteria found in soil.
- ▶ They are characterized by their tough, leathery, frequently pigmented colonies and their filamentous growth.
- ▶ *Streptomyces* are chemoheteroorganotrophs, growing best at 25°C and pH 8–9.
- ▶ They use complex organic materials as carbon and energy sources and are involved in the breakdown of these products in the soil.
- ▶ This degradative ability makes these bacteria pivotal in the production of fertile soil for agriculture.

- ▶ They also give soil its characteristic **smell** by the production of volatile low molecular weight compounds called geosmins.
- ▶ *Streptomyces* are also of medical and industrial importance because they synthesize **antibiotics**.
- ▶ Those antibiotics help the organism compete with other organisms in the nutrient-depleted environment of the soil by reducing competition.
- ▶ Over 50 different antibiotics have been isolated from *Streptomyces* species, including **streptomycin**, **neomycin**, **chloramphenicol** and **tetracyclines**.



# Identification and Examination of soil micro-organisms


## Materials

- ▶ 1–2 media plates
- ▶ wire loop
- ▶ Bunsen burner

***Streptomyces* species are white and colorful chalky looking colonies.**



# 1 – Macroscopic examination



- ▶ Examine the ISP2 plates and look for typical *Streptomyces* colonies.
  - ▶ They are small, opaque, compact, frequently pigmented (brown, yellow, pink, etc.), often leathery, and appear dry and dull looking.
  - ▶ Typically, a depression in the agar surface will be observed around the colony.
  - ▶ Avoid molds. They usually form much softer, fuzzy colonies if present.
- 

- ▶ Good *Streptomyces* candidates will be difficult to remove from the agar with the inoculating needle or loop and upon observation under the microscope will reveal a multitude of spores with a few filamentous cells.




## 2– Staining of *Sterptomyces* by Gram stain

- ▶ Using a flamed loop, scrape part of a **colony** onto a drop of **water** on a **slide**. Let dry.
- ▶ Heat fix by gently passing the slide through the flame 3 or 4 times.
- ▶ Flood slide with **crystal violet** for 15 secs, wash with water.
- ▶ Flood with **iodine** for 1 min, wash.
- ▶ Add 2 drops of **95% ethanol** for 5 seconds, wash.
- ▶ Flood with **safranin** for 5 min, wash.
- ▶ Dry the slide, add 1 drop of cedar wood oil.
- ▶ View using oil immersion lens (100X).

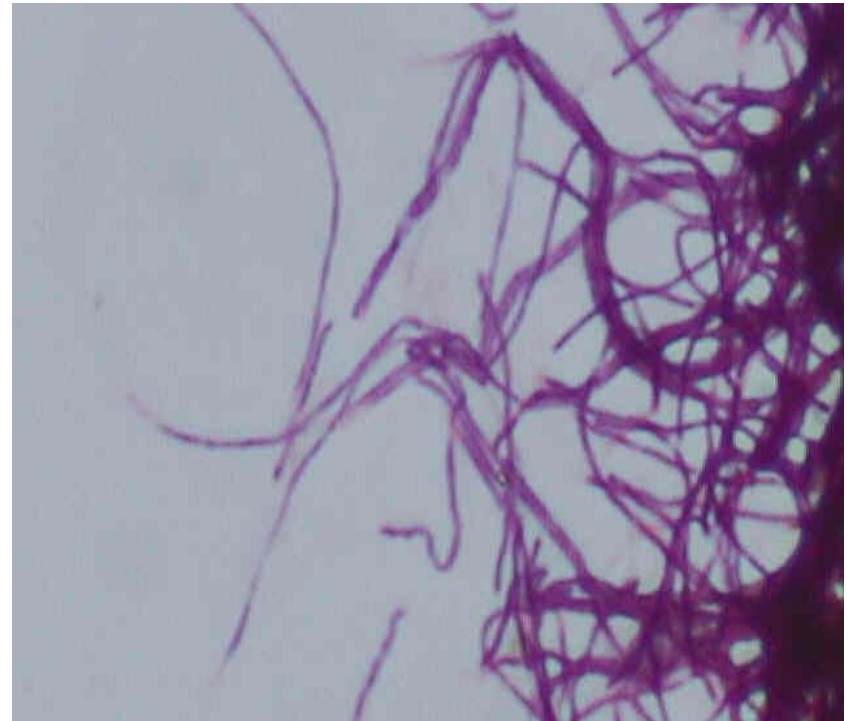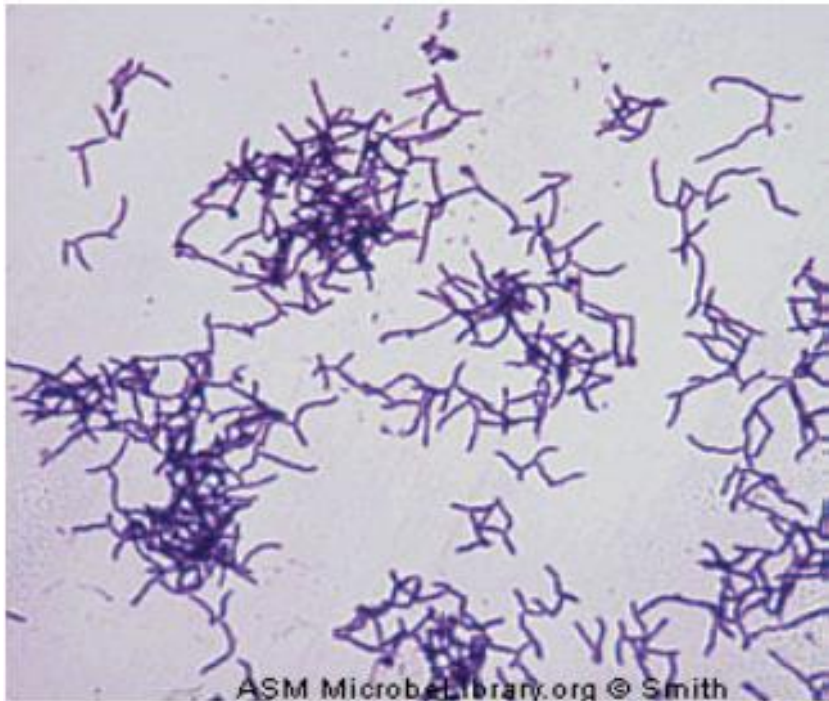


# Scheme for Streptomyces stained by Gram stain

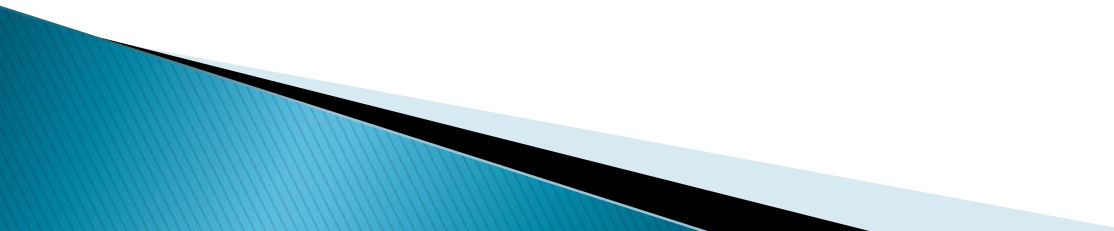
- ▶ **Type of stain**: Differential stain (Gram stain).
  - ▶ **Name of dye**: Crystal violet (1<sup>st</sup> stain), Iodine (mordant), Ethanol 95% (decolorizer) and Safranin (counter stain)
  - ▶ **Color**: Gram positive stained with deep violet
  - ▶ **Culture**: single.
  - ▶ **Shape**: Threads.
- 

## ❑ Comment


- ▶ In Gram **positive** bacteria, cell wall contain **low** lipid content.
- ▶ 95% ethanol dissolves lipids , leaving **small pores not allowing** the escape of (Iodine–crystal violet complex).
- ▶ So the bacteria retain the deep violet colour.
  
- ▶ In Gram **negative** bacteria, cell wall contain **high** lipid content.
- ▶ 95% ethanol dissolves lipids , leaving **large pores allowing** the escape of (Iodine–crystal violet complex).
- ▶ So when saffranin is added ,the bacteria stained with red colour.

# 3– Streaking for isolation

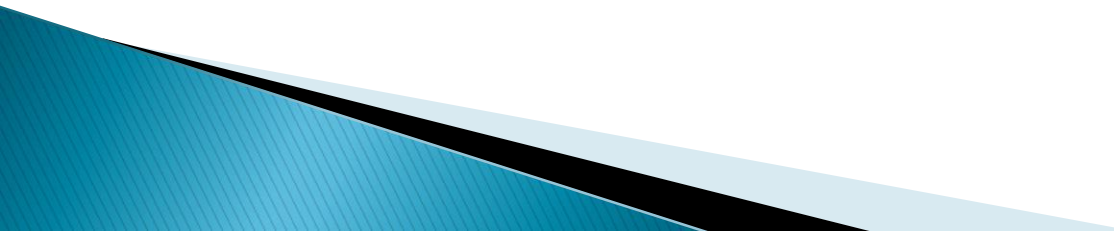- ▶ With an inoculating loop streak *Streptomyces* colony on the transfer media for isolation of pure colonies.
  - ▶ Incubate for 3–5 days at room temp (at dark)
- 

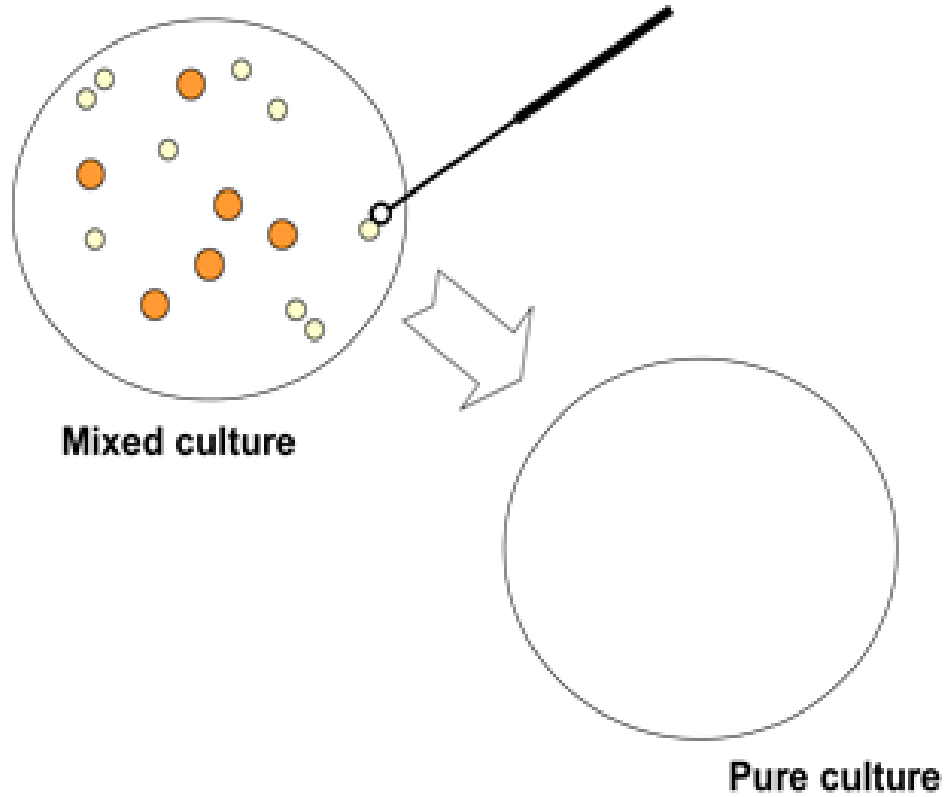



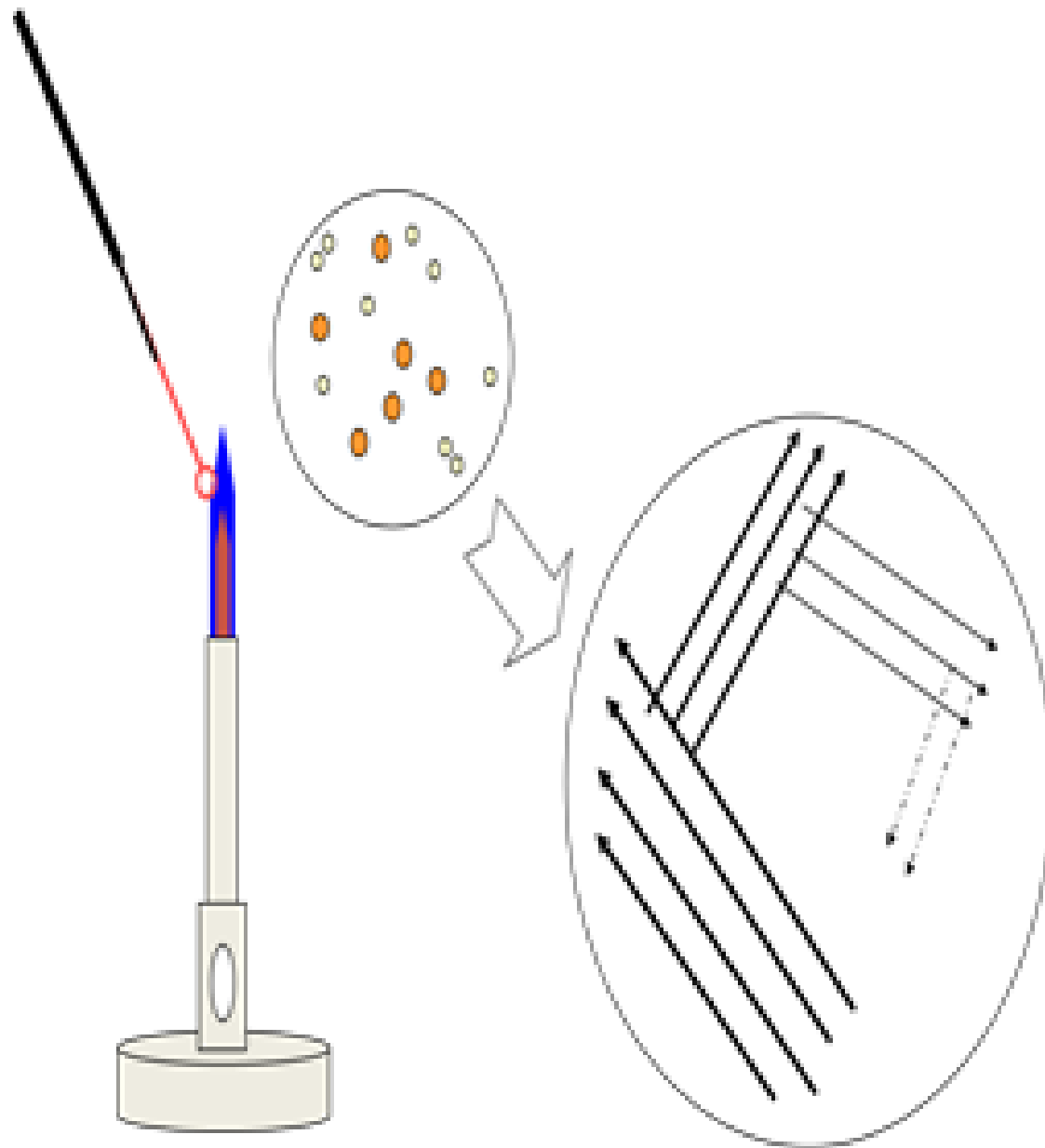

# After incubation



*Thank  
you*

